

# neptune

NEW



- Robust, cost effective, easy to re-issue
- Easy to adjust and set up for use
- Vast range of adjustment for growth




# There are 2 seat sizes available to cover a vast age range of 18 months to 10 years approx.



Adjustable Headrest



Flexi or Lateral Supports

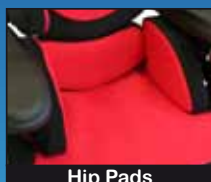

Hip Pads



Tray




Footplate & Leg Rest



Footplate Pad



Sandals or Ankle Supports

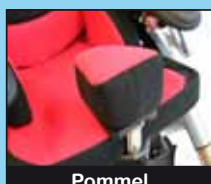

Pommel



Contoured Seat Cushion



Positioning Aids

A user friendly design ensures that the Neptune can be easily set up and adjusted during use.

Referencing the EN 1729 standard for height adjustment in classroom furniture, the Neptune has an excellent height range; specially designed for



nursery, school and home environments.

The Neptunes bio-mechanical reclining back rotates around the hip joint, ensuring the trunk supports remain in the optimum position.

The robust, yet lightweight design ensures effortless manoeuvrability and compact transportation.

The Neptune also has removable upholstery which is easy to clean or replace so the



product can be re-issued without fuss.

A strong, adjustable one piece footplate can be removed or flipped up for standing transfers.

## Stain Resistant Fabrics



red/black



blue/black



orange/black



pink/black

NEPTUNE (mm)	SIZE 1		SIZE 2	
Measurements	Min	Max	Min	Max
Seat Depth	200	325	325	450
Seat to Top Back	325	450	425	575
Width at Hips	165	325	265	425
Width at Chest	150	270	225	345
Seat to Footplate	130	290	260	420
Weight of Seating Unit	10kg		13kg	
Indoor Base Frame				
Seat to Floor: Min 240 - Max 540				
Tilt Angle: 25° / Prone Angle: 15°				

